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
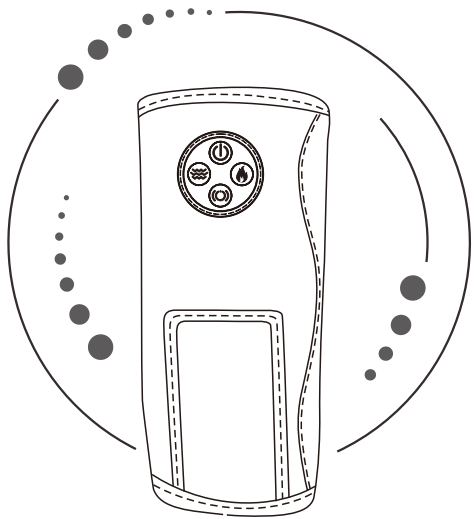
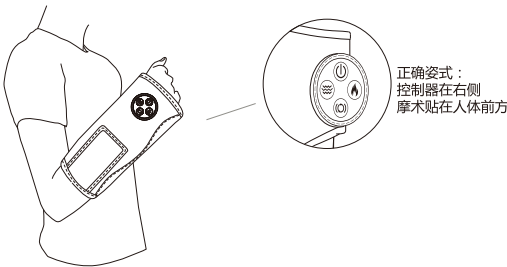
ARM (CALF) MASSAGER

air pressure / heating / viberation /
Red light roasting / Timing massage /

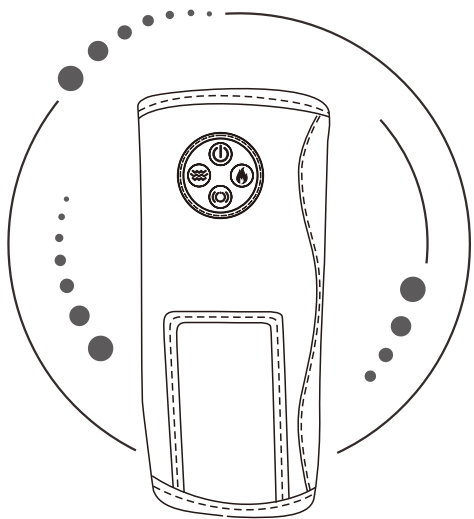operation instruction and user manual


PRODUCT INSTRUCTION AND FUNCTION:

- ①The 5V2A USB conn ector is plugged in with a power bank or a 5V2Aadapteror DC 5V2Ain car
- ②easily carry and enjoy massage.
- ③product surface is good looking and fashion, leading trend. .
- ④Red wavelength 650nm effectively alleviatediscomfort, and both traditional and modern science technology multi-functional massage
- ⑤Maxium rated power 8W and working long time.


USE PRODUCT METHOD

1、pre-use preparation

- (1)take it out

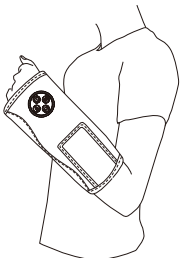
- (2)connect battery bank or adaptor,press the power key 2 sec

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
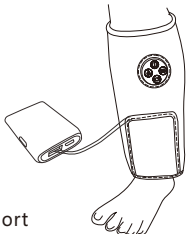
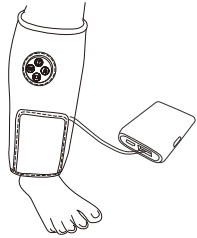
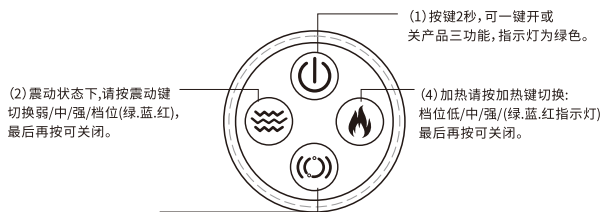
2、operation instruction (below picture)

- (1) open magic sticker,bundle it to comfortable positlon

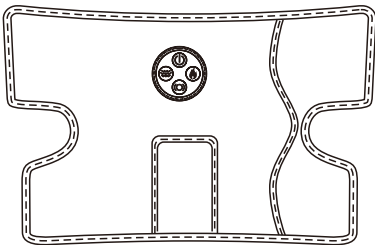- (2) use magic sticker to adjust right position

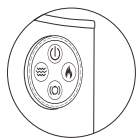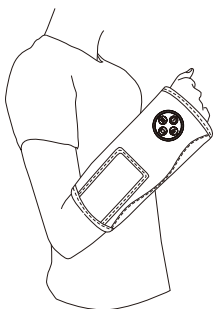right position:
controller on top half,
battery bank pocket on lower half.

- (3)Press power key 2 sec. All function working.green light open
- (4) Red light heating state:low/middle/ High /off(light color:Green/ blue/ red)
- (5)Stereoscopic wave airbag state:I / II /III mode/off(light color:Green/ blue/ red)
- (6)Vibration mode state:weak/pulse/steady/off (light color:Green/ blue/ red)
- (7)automatic off when continuous working 15 minutes.

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
3、controller instruction

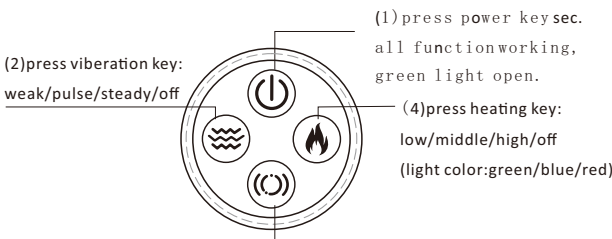- (3)press air pressure key:
I/II/III mode/off(Simulated artificial massage)

- (4)A heating chip and vibrating motor are built into the bottom of the foot

4、precautions!

- (1)Be careful to use the highest heating for long periods oftime to avoid burning skin.
- (2)In order to prolong the service life of the product and avoid long-term unlimited use, it is recommended to use the product no more than three times a day for 15 minutes each time.
- (3)prohibited use product when sleeping
- (4)can wipe surface of product,prohibited into waterto clean

5.disable crowds:

- (1) skin with thermal allergies
- (2)unable take care of themselves
- (3)doctor warning